




J DENTAL CARE

just smile

JD EVOLUTION ^{plus}

Simple insertion, natural stability

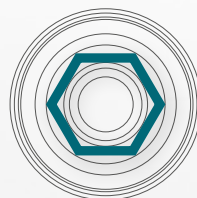

hexagonal connection




Made in Italy

SIMPLE INSERTION, NATURAL STABILITY


features & advantages

- Internal hex connection with lead-in bevel
- Anatomic coronal collar
- Expanding tapered inner body
- Self-cutting capacity in both directions
- Apical cutting blades
- The connection is universally recognized by all the main digital libraries
- Reduces stress concentration on the crestal bone and maximizes bone volume around the implant neck for stable hard and soft tissues
- Makes it possible to achieve high primary stability

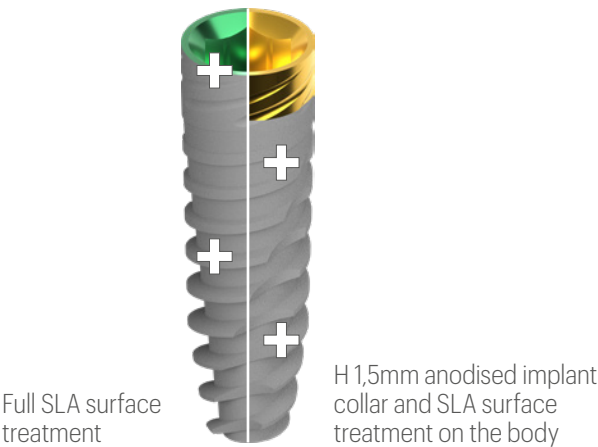

SINGLE PROSTHETIC CONNECTION

Single prosthetic connection for all diameters
Reduced inventory and enhanced working flexibility

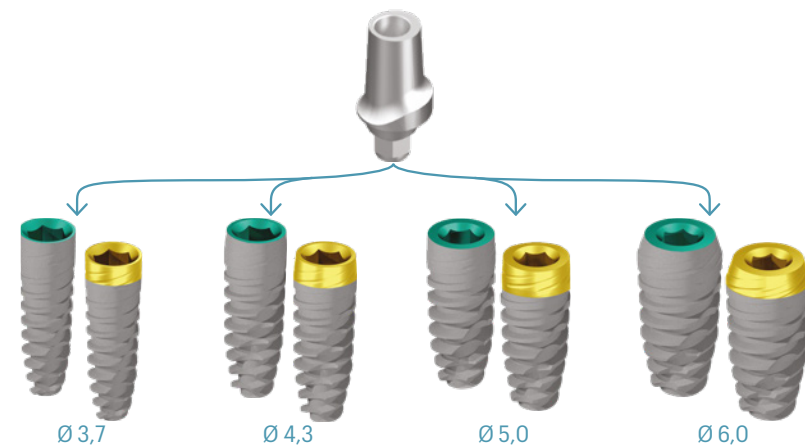

PRODUCT SPECIFICATIONS


Note: All measurements in mm

The JDEvolution Plus dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	LENGTH	TIP DIAMETER	PLATFORM DIAMETER	ABUTMENT INTERFACE
Ø 3.7	8	2.8	3.7	3.4
	10, 11.5, 13, 15, 18	2.4		
Ø 4.0	6	3.2	4.0	3.4
	8	2.8		
	10	2.4		
	11.5, 13	2.4		
Ø 4.3	6, 8	3.6	4.0	3.4
	10, 11.5, 13, 15, 18	2.8		
Ø 5.0	6	3.6	4.8	3.4
	8, 10, 11.5, 13, 15	3.2		
Ø 6.0	6	4.4	5.0	3.4
	8, 10, 11.5, 13, 15	4.0		

The JDEvolution Plus T dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	LENGTH	TIP DIAMETER	PLATFORM DIAMETER	ABUTMENT INTERFACE
Ø 4.0	6, 8, 10, 11.5, 13, 15, 18	2.4	3.8	3.4

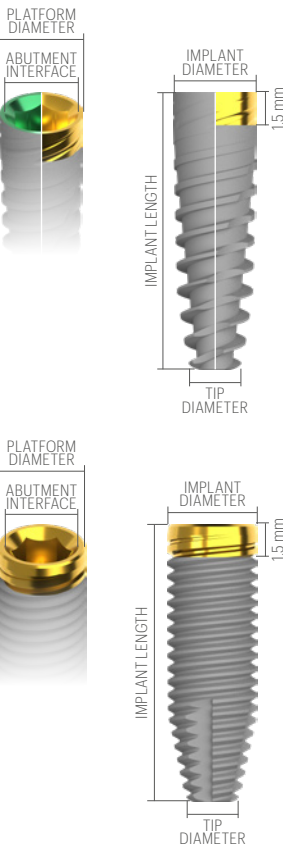

PRODUCT CATALOGUE

Note: all measurements of diameters and lengths are in mm

Implants:

Ø 3,7

EV37080:	JDEvolution Plus Ø 3.7 L 8
EV37100:	JDEvolution Plus Ø 3.7 L 10
EV37115:	JDEvolution Plus Ø 3.7 L 11.5
EV37130:	JDEvolution Plus Ø 3.7 L 13
EV37150:	JDEvolution Plus Ø 3.7 L 15
EV37180:	JDEvolution Plus Ø 3.7 L 18

Ø 4,0

EV40060:	JDEvolution Plus Ø 4.0 L 6
EV40080:	JDEvolution Plus Ø 4.0 L 8
EV40100:	JDEvolution Plus Ø 4.0 L 10
EV40115:	JDEvolution Plus Ø 4.0 L 11.5
EV40130:	JDEvolution Plus Ø 4.0 L 13

Ø 4,3

EV43060:	JDEvolution Plus Ø 4.3 L 6
EV43080:	JDEvolution Plus Ø 4.3 L 8
EV43100:	JDEvolution Plus Ø 4.3 L 10
EV43115:	JDEvolution Plus Ø 4.3 L 11.5
EV43130:	JDEvolution Plus Ø 4.3 L 13
EV43150:	JDEvolution Plus Ø 4.3 L 15
EV43180:	JDEvolution Plus Ø 4.3 L 18

Ø 5,0

EV50060:	JDEvolution Plus Ø 5.0 L 6
EV50080:	JDEvolution Plus Ø 5.0 L 8
EV50100:	JDEvolution Plus Ø 5.0 L 10
EV50115:	JDEvolution Plus Ø 5.0 L 11.5
EV50130:	JDEvolution Plus Ø 5.0 L 13
EV50150:	JDEvolution Plus Ø 5.0 L 15
EV50180:	JDEvolution Plus Ø 5.0 L 18

Ø 6,0

EV60060:	JDEvolution Plus Ø 6.0 L 6
EV60080:	JDEvolution Plus Ø 6.0 L 8
EV60100:	JDEvolution Plus Ø 6.0 L 10
EV60115:	JDEvolution Plus Ø 6.0 L 11.5
EV60130:	JDEvolution Plus Ø 6.0 L 13
EV60150:	JDEvolution Plus Ø 6.0 L 15



Anodised implant collar:

Ø 3,7


EVA37080:	JDEvolution Plus Anodized Collar Ø 3.7 L 8
EVA37100:	JDEvolution Plus Anodized Collar Ø 3.7 L 10
EVA37115:	JDEvolution Plus Anodized Collar Ø 3.7 L 11.5
EVA37130:	JDEvolution Plus Anodized Collar Ø 3.7 L 13
EVA37150:	JDEvolution Plus Anodized Collar Ø 3.7 L 15
EVA37180:	JDEvolution Plus Anodized Collar Ø 3.7 L 18

Ø 4,0

EVA40150:	JDEvolution Plus Anodized Collar Ø 4.0 L 15
EVA40180:	JDEvolution Plus Anodized Collar Ø 4.0 L 18

Ø 4,3

EVA43060:	JDEvolution Plus Anodized Collar Ø 4.3 L 6
EVA43080:	JDEvolution Plus Anodized Collar Ø 4.3 L 8
EVA43100:	JDEvolution Plus Anodized Collar Ø 4.3 L 10
EVA43115:	JDEvolution Plus Anodized Collar Ø 4.3 L 11.5
EVA43130:	JDEvolution Plus Anodized Collar Ø 4.3 L 13
EVA43150:	JDEvolution Plus Anodized Collar Ø 4.3 L 15
EVA43180:	JDEvolution Plus Anodized Collar Ø 4.3 L 18

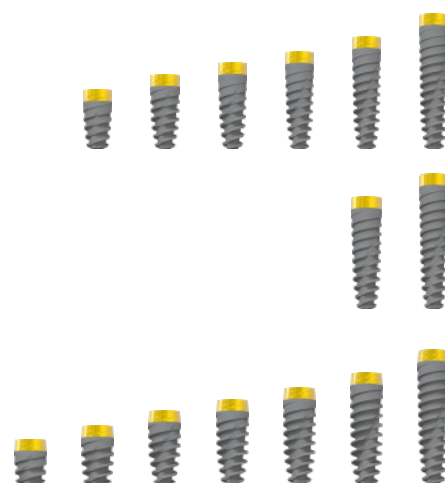

Ø 5,0

EVA50060:	JDEvolution Plus Anodized Collar Ø 5.0 L 6
EVA50080:	JDEvolution Plus Anodized Collar Ø 5.0 L 8
EVA50100:	JDEvolution Plus Anodized Collar Ø 5.0 L 10
EVA50115:	JDEvolution Plus Anodized Collar Ø 5.0 L 11.5
EVA50130:	JDEvolution Plus Anodized Collar Ø 5.0 L 13
EVA50150:	JDEvolution Plus Anodized Collar Ø 5.0 L 15
EVA50180:	JDEvolution Plus Anodized Collar Ø 5.0 L 18

Ø 6,0

EVA60060:	JDEvolution Plus Anodized Collar Ø 6.0 L 6
EVA60080:	JDEvolution Plus Anodized Collar Ø 6.0 L 8
EVA60100:	JDEvolution Plus Anodized Collar Ø 6.0 L 10
EVA60115:	JDEvolution Plus Anodized Collar Ø 6.0 L 11.5
EVA60130:	JDEvolution Plus Anodized Collar Ø 6.0 L 13
EVA60150:	JDEvolution Plus Anodized Collar Ø 6.0 L 15



JDEvolution Plus T:

EVT40060:	JDEvolution Plus T Anodized Collar Ø 4.0 L 6
EVT40080:	JDEvolution Plus T Anodized Collar Ø 4.0 L 8
EVT40100:	JDEvolution Plus T Anodized Collar Ø 4.0 L 10
EVT40115:	JDEvolution Plus T Anodized Collar Ø 4.0 L 11.5
EVT40130:	JDEvolution Plus T Anodized Collar Ø 4.0 L 13
EVT40150:	JDEvolution Plus T Anodized Collar Ø 4.0 L 15
EVT40180:	JDEvolution Plus T Anodized Collar Ø 4.0 L 18



Cover screw:

EVCS:	Cover Screw JDEvolution Plus
EVCS6:	Cover Screw Ø 6.0 JDEvolution Plus



Healing Abutment:


EVHA41:	Healing Abutment Ø 4.0 H 1.0 JDEvolution Plus
EVHA42:	Healing Abutment Ø 4.0 H 2.0 JDEvolution Plus
EVHA43:	Healing Abutment Ø 4.0 H 3.0 JDEvolution Plus
EVHA45:	Healing Abutment Ø 4.0 H 5.0 JDEvolution Plus
EVHA47:	Healing Abutment Ø 4.0 H 7.0 JDEvolution Plus
EVHA53:	Healing Abutment Ø 5.0 H 3.0 JDEvolution Plus
EVHA55:	Healing Abutment Ø 5.0 H 5.0 JDEvolution Plus
EVHA57:	Healing Abutment Ø 5.0 H 7.0 JDEvolution Plus
EVHA59:	Healing Abutment Ø 5.0 H 9.0 JDEvolution Plus
EVHA63:	Healing Abutment Ø 6.0 H 3.0 JDEvolution Plus
EVHA65:	Healing Abutment Ø 6.0 H 5.0 JDEvolution Plus
EVHA67:	Healing Abutment Ø 6.0 H 7.0 JDEvolution Plus
EVHA753:	Healing Abutment Ø 7.5 H 3.0 JDEvolution Plus
EVHA755:	Healing Abutment Ø 7.5 H 5.0 JDEvolution Plus
EVHA54TS:	Healing Abutment Tissue Ø 5.0 H 4.0 JDEvolution Plus
EVHA56TS:	Healing Abutment Tissue Ø 5.0 H 6.0 JDEvolution Plus



Impression coping:

OPEN TRAY

EVICOT4C:	Impression Coping Open Tray Ø 4.0 JDEvolution Plus
EVICOT5C:	Impression Coping Open Tray Ø 5.0 JDEvolution Plus
EVICOT6C:	Impression Coping Open Tray Ø 6.0 JDEvolution Plus
EVICOTNEC:	Impression Coping Open Tray Non Engaging JDEvolution Plus
EVICOTS:	Impression Coping Open Tray Screw JDEvolution Plus
EVICOT5TSC:	Impression Coping Tissue Open Tray Ø 5.0 JDEvolution Plus

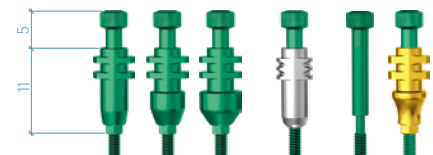

CLOSED TRAY


EVICCT4C:	Impression Coping Closed Tray Ø 4.0 JDEvolution Plus
EVICCT5C:	Impression Coping Closed Tray Ø 5.0 JDEvolution Plus
EVICCT6C:	Impression Coping Closed Tray Ø 6.0 JDEvolution Plus
EVICCTS:	Impression Coping Closed Tray Screw JDEvolution Plus
EVICCT1C:	Impression Coping Closed Tray Direct Press Fit JDEvolution Plus



Temporary Abutment:

EVTAE:	Temporary Abutment Engaging JDEvolution Plus
EVTAEPC:	Temporary Abutment Engaging Peek JDEvolution Plus
EVTANEC:	Temporary Abutment Non Engaging JDEvolution Plus
EVTATSC:	Temporary Abutment Tissue Engaging JDEvolution Plus
EVTANETSC:	Temporary Abutment Tissue Non Engaging JDEvolution Plus

Torque recommended 30 Ncm

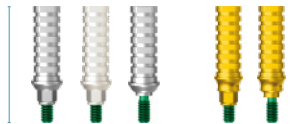

Straight Abutment:

EVNSA4510C:	Straight Abutment Ø 4.5 H 1.5 JDEvolution Plus
EVNSA4520C:	Straight Abutment Ø 4.5 H 2.0 JDEvolution Plus
EVNSA4540C:	Straight Abutment Ø 4.5 H 4.0 JDEvolution Plus
EVNSA5020C:	Straight Abutment Ø 5.0 H 2.0 JDEvolution Plus
EVNSA5040C:	Straight Abutment Ø 5.0 H 4.0 JDEvolution Plus
EVNSA6020C:	Straight Abutment Ø 6.0 H 2.0 JDEvolution Plus
EVNSA6040C:	Straight Abutment Ø 6.0 H 4.0 JDEvolution Plus
EVNSA5020TSC:	Straight Abutment Tissue Ø 5.0 H 2.0 JDEvolution Plus
EVNSA5040TSC:	Straight Abutment Tissue Ø 5.0 H 4.0 JDEvolution Plus

Torque recommended 30 Ncm



Torque recommended 30 Ncm



GP Abutment:

EV15GPAEC:	GP Abutment 15° JDEvolution Plus
EV25GPAEC:	GP Abutment 25° JDEvolution Plus
EVGPA40EC:	GP Abutment Engaging Ø 4.0 JDEvolution Plus
EVGPA40NEC:	GP Abutment Non Engaging Ø 4.0 JDEvolution Plus
EVGPAEC:	GP Abutment Engaging Ø 5.0 JDEvolution Plus
EVGPA50NEC:	GP Abutment Non Engaging Ø 5.0 JDEvolution Plus
EVGPAEPC:	GP Abutment Engaging Peek Ø 5.0 JDEvolution Plus
EVGPA60EC:	GP Abutment Engaging Ø 6.0 JDEvolution Plus

Torque recommended 30 Ncm




Torque recommended 30 Ncm



Anatomic Abutment:

EVEA15:	Anatomic Abutment H 1.5 JDEvolution Plus
EVEA30:	Anatomic Abutment H 3.0 JDEvolution Plus
EV15A15:	15° Angulated Anatomic Abutment Ø 5.0 H 1.5 JDEvolution Plus
EV15A15N:	15° Angulated Anatomic Abutment Ø 4.5 H 1.5 JDEvolution Plus
EV15A30:	15° Angulated Anatomic Abutment Ø 5.0 H 3.0 JDEvolution Plus
EV15A30N:	15° Angulated Anatomic Abutment Ø 4.5 H 3.0 JDEvolution Plus
EV25A15:	25° Angulated Anatomic Abutment Ø 5.0 H 1.5 JDEvolution Plus
EV25A30:	25° Angulated Anatomic Abutment Ø 5.0 H 3.0 JDEvolution Plus

Torque recommended 30 Ncm

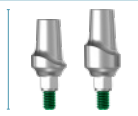

Torque recommended 30 Ncm



Rapid Abutment:

EVRA15C:	Rapid Abutment H 1.5 JDEvolution Plus
EVRA30C:	Rapid Abutment H 3.0 JDEvolution Plus
EVRA15D:	Rapid Abutment w/o Healing Cap and Coping H 1.5 JDEvolution Plus
EVRA30D:	Rapid Abutment w/o Healing Cap and Coping H 3.0 JDEvolution Plus
EVRAHC:	Rapid Abutment Healing Cap
EVRAIC:	Rapid Abutment Impression Coping
EVRAAR:	Rapid Abutment Replica

Torque recommended 30 Ncm




Torque recommended 30 Ncm



Wax Up Abutment:

EVWAEC:	Wax-Up Abutment Engaging JDEvolution Plus
EVWANEC:	Wax-Up Abutment Non Engaging JDEvolution Plus

Torque recommended 30 Ncm

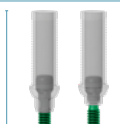

Cobalt-Chrome Abutment:

EVCCCEC:	Chrome-Cobalt Abutment Engaging JDEvolution Plus
EVCCNEC:	Chrome-Cobalt Abutment Non Engaging JDEvolution Plus
EVCAACCNEC:	Cobalt-Chrome Conical Abutment Non Engaging JDEvolution Plus

Torque recommended 30 Ncm



15 Ncm

Gold Abutment:

EVGAEC:	Gold Abutment Engaging JDEvolution Plus
EVGANEC:	Gold Abutment Non Engaging JDEvolution Plus

Torque recommended 30 Ncm



Conical Abutment:

EVNCA10C:	Conical Abutment Straight H 1.0 JDEvolution Plus
EVNCA15C:	Conical Abutment Straight H 1.5 JDEvolution Plus
EVNCA20C:	Conical Abutment Straight H 2.0 JDEvolution Plus
EVNCA30C:	Conical Abutment Straight H 3.0 JDEvolution Plus
EVNCA40C:	Conical Abutment Straight H 4.0 JDEvolution Plus
EVNCA70C:	Conical Abutment Straight H 7.0 JDEvolution Plus
EVNCA90C:	Conical Abutment Straight H 9.0 JDEvolution Plus
EVNCA15TSC:	Conical Abutment Tissue Straight H 1.5 JDEvolution Plus
EVNCA30TSC:	Conical Abutment Tissue Straight H 3.0 JDEvolution Plus
EVNCA40TSC:	Conical Abutment Tissue Straight H 4.0 JDEvolution Plus
EVNCA50TSC:	Conical Abutment Tissue Straight H 5.0 JDEvolution Plus
EVCA1725C:	Conical Abutment Angulated 17° H 2.5 JDEvolution Plus
EVCA1735C:	Conical Abutment Angulated 17° H 3.5 JDEvolution Plus
EVCA3025C:	Conical Abutment Angulated 30° H 2.5 JDEvolution Plus
EVCA3035C:	Conical Abutment Angulated 30° H 3.5 JDEvolution Plus
EVCA1725TSC:	Conical Abutment Tissue Angulated 17° H 2.5 JDEvolution Plus
EVCA1735TSC:	Conical Abutment Tissue Angulated 17° H 3.5 JDEvolution Plus
EVA275C:	Conical Abutment Tissue Angulated 17° H 5.0 JDEvolution Plus
EVCA3025TSC:	Conical Abutment Tissue Angulated 30° H 2.5 JDEvolution Plus
EVCA3035TSC:	Conical Abutment Tissue Angulated 30° H 3.5 JDEvolution Plus
EVCA3050TSC:	Conical Abutment Tissue Angulated 30° H 5.0 JDEvolution Plus
EVA137C:	Conical Abutment Straight Engaging H 1.5 JDEvolution Plus
EVA138C:	Conical Abutment Straight Engaging H 3.0 JDEvolution Plus
EVS105:	Screw for Conical Abutment Straight Engaging H 1.5 JDEvolution Plus
EVS106:	Screw for Conical Abutment Straight Engaging H 3.0 JDEvolution Plus
EVCASA:	Screw for Angulated Conical Abutment JDEvolution Plus
EVCAPS:	Prosthetic Screw Conical Abutment JDEvolution Plus
EVCAPSL:	Prosthetic Screw Conical Abutment Long JDEvolution Plus
EVCAPSA	Prosthetic Screw Angulated for Conical Abutment JDEvolution *
EVCAPSN:	Prosthetic Screw Conical Abutment Direct N
EVCAPSN A	Prosthetic Screw for Angulated Channel Conical Abutment Direct N *
EVCAICOTEC:	Conical Abutment Impression Coping Engaging Open Tray JDEvolution Plus
EVCAICOTC:	Conical Abutment Impression Coping Open Tray JDEvolution Plus
EVCAICOTLC:	Conical Abutment Impression Coping Open Tray Long JDEvolution Plus
EVCAICOT02:	Screw for Conical Abutment Impression Coping Open Tray JDEvolution Plus
EVCAICOT04:	Screw for Conical Abutment Impression Coping Open Tray Long JDEvolution Plus
EVCAICCTC	Conical Abutment Impression Coping Closed Tray
EVCAHC:	Healing Cap for Conical Abutment JDEvolution Plus
EVCAHCL:	Conical Abutment Healing Cap H 6.0 JDEvolution Plus
EVCAHC9:	Conical Abutment Healing Cap H 9.0 JDEvolution Plus
EVCAHCB:	Conical Abutment Healing Cap Peek JDEvolution Plus
EVCAGPAEC:	GP Abutment Engaging for Conical Abutment JDEvolution Plus
EVCAGPANEC:	GP Abutment Non Engaging for Conical Abutment JDEvolution Plus
EVCATANEWSC:	Temporary Abutment Non Engaging Conical Abutment Smooth for Welding JDEvolution Plus
EVCATANEC:	Temporary Abutment Non Engaging Conical Abutment JDEvolution Plus
EVCATANEWNC:	Temporary Abutment Non Engaging Conical Abutment for Welding JDEvolution Plus
EVCAWANEC:	Wax-Up Abutment Non Engaging for Conical Abutment JDEvolution Plus
EVCAAR	Conical Abutment Replica

* for this screw use Angulated Screw Driver for JDTorque




Torque recommended 30 Ncm



Torque recommended 30 Ncm



Torque recommended 30 Ncm

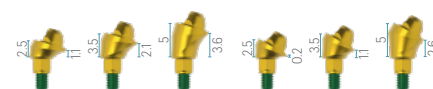


Torque recommended 30 Ncm



Torque recommended 30 Ncm



Torque recommended 30 Ncm

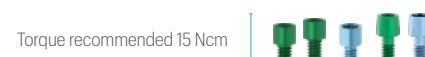

Torque recommended 15 Ncm



Torque recommended 15 Ncm




Torque recommended 15 Ncm

Conical Abutments Mini:

EVNCA10MIC:	Conical Straight Abutment MINI H 1.0 JDEvolution Plus
EVNCA15MIC:	Conical Straight Abutment MINI H 1.5 JDEvolution Plus
EVNCA30MIC:	Conical Straight Abutment MINI H 3.0 JDEvolution Plus
EVNCA40MIC:	Conical Straight Abutment MINI H 4.0 JDEvolution Plus
EVNCA50MIC:	Conical Straight Abutment MINI H 5.0 JDEvolution Plus
EVCAHCMIC:	Conical Abutment MINI Healing Cap JDEvolution Plus
EVCAHCLMIC:	Conical Abutment MINI Healing Cap H 6.0 JDEvolution Plus
EVCAITMIC:	Interface CAD CAM Non Engaging H 7.0 for Conical Abutment MINI JDEvolution Plus
EVCAITMIEC:	Interface CAD CAM Engaging H 7.0 for Conical Abutment MINI JDEvolution Plus
EVCATIBMIC:	Interface CAD CAM Non Engaging H 5.0 for Conical Abutment MINI JDEvolution Plus
EVCATIBMIEC:	Interface CAD CAM Engaging H 5.0 for Conical Abutment MINI JDEvolution Plus



EVCATANEMIC:	Temporary Abutment Non Engaging Conical Abutment MINI JDEvolution Plus
EVCAICOTMIC:	Conical Abutment MINI Impression Coping Open Tray JDEvolution Plus
EVCAICOTMIEC:	Conical Abutment MINI Impression Coping Engaging Open Tray JDEvolution Plus
EVCAARMI	Conical Abutment MINI Replica
EVCAARMIC	Conical Abutment MINI Replica CAD CAM
EVCASBCEMIC:	ScanBody Conical Abutment MINI JDEvolution Plus
EVCASBCEMIEC:	ScanBody Engaging Conical Abutment MINI JDEvolution Plus
EVCASBCESMIC:	On Top ScanBody Conical Abutment MINI

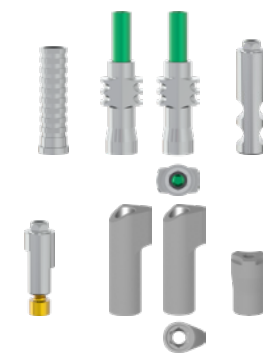

Prosthetic screw:

EVS:	Prosthetic Screw JDEvolution Plus
EVSA:	Prosthetic Screw Angulated JDEvolution Plus *
* for this screw use Angulated Screw Driver for JD Torque	


Torque recommended 30 Ncm



Interfaces:

EVINTEC:	Interface Engaging JDEvolution Plus
EVINTNEC:	Interface Non Engaging JDEvolution Plus
70287C:	Interface Non Engaging for Conical Abutment JDEvolution Plus


Torque recommended 30 Ncm

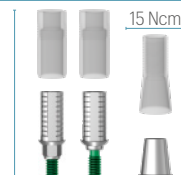

AS Abutment:


EVASCEC:	AS Abutment Engaging JDEvolution Plus
EVASCNEC:	AS Abutment Non Engaging JDEvolution Plus

Torque recommended 30 Ncm

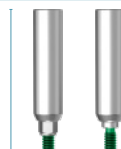

TiBase, JDScanBody and Implant Replica CAD CAM:

EVSBCEC:	ScanBody Engaging JDEvolution Plus
EVCASBCEC:	ScanBody Conical Abutment JDEvolution Plus
EVCASBCEEC:	ScanBody Engaging Conical Abutment JDEvolution Plus
EVCASBCECSC:	On Top ScanBody Conical Abutment JDEvolution Plus
JDSCGIOSE	ScanBody Full Arch Tracer Hexagonal
JDSCGIOS	ScanBody Full Arch Tracer
JDSCGIOSR	ScanBody Rectangular Full Arch Tracer

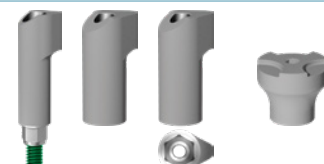

EVHSB53:	On Top ScanBody Engaging Ø 5.0 H 3.0 JDEvolution Plus
EVHSB55:	On Top ScanBody Engaging Ø 5.0 H 5.0 JDEvolution Plus
EVHSB57:	On Top ScanBody Engaging Ø 5.0 H 7.0 JDEvolution Plus
EVSCPC:	ScanPost JDEvolution Plus



EVTIBC:	TiBase Engaging Ø4.25 H 0.5 C 4.7 JDEvolution Plus
EVTIBNEC:	TiBase Non Engaging Ø4.25 H 0.5 C 4.7 JDEvolution Plus
EVTIB15C:	TiBase Engaging Ø4.25 H 1.5 C 4.7 JDEvolution Plus
EVTIB15NEC:	TiBase Non Engaging Ø4.25 H 1.5 C 4.7 JDEvolution Plus
EVTIB30C:	TiBase Engaging Ø4.25 H 3.0 C 4.7 JDEvolution Plus
EVTIB30NEC:	TiBase Non Engaging Ø4.25 H 3.0 C 4.7 JDEvolution Plus

Torque recommended 30 Ncm



EVTIBNC:	TiBase Long Engaging Ø4.4 H 0.5 C 7.7 JDEvolution Plus
EVTIBNENC:	TiBase Long Non Engaging Ø4.4 H 0.5 C 7.7 JDEvolution Plus
EVTIB10NC:	TiBase New Engaging Ø4.4 H 1.0 C 7.7 JDEvolution Plus
EVTIB10NENC:	TiBase New Non Engaging Ø4.4 H 1.0 C 7.7 JDEvolution Plus
EVTIB15NC:	TiBase Long Engaging Ø4.4 H 1.5 C 7.7 JDEvolution Plus
EVTIB15NENC:	TiBase Long Non Engaging Ø4.4 H 1.5 C 7.7 JDEvolution Plus
EVTIB20NC:	TiBase New Engaging Ø4.4 H 2.0 C 7.7 JDEvolution Plus
EVTIB20NENC:	TiBase New Non Engaging Ø4.4 H 2.0 C 7.7 JDEvolution Plus
EVTIB30NC:	TiBase Long Engaging Ø4.4 H 3.0 C 7.7 JDEvolution Plus
EVTIB30NENC:	TiBase Long Non Engaging Ø4.4 H 3.0 C 7.7 JDEvolution Plus

Torque recommended 30 Ncm



EVA150C:	TiBase Long Engaging Ø4.4 H 0.5 C 7.7 for Angulated Channel JDEvolution Plus *
EVA151C:	TiBase New Non Engaging Ø4.4 H 0.5 C 7.7 for Angulated Channel JDEvolution Plus *
EVA152C:	TiBase Long Engaging Ø4.4 H 1.5 C 7.7 for Angulated Channel JDEvolution Plus *
EVA153C:	TiBase New Non Engaging Ø4.4 H 1.5 C 7.7 for Angulated Channel JDEvolution Plus *
EVA154C:	TiBase Long Engaging Ø4.4 H 3.0 C 7.7 for Angulated Channel JDEvolution Plus *
EVA155C:	TiBase New Non Engaging Ø4.4 H 3.0 C 7.7 for Angulated Channel JDEvolution Plus *

Torque recommended 30 Ncm



EVA133C:	Tibase Interface Engaging H 6.0 for Angulated Channel JDEvolution Plus **	Torque recommended 30 Ncm	
EVA134C:	Tibase Interface Non Engaging H 6.0 for Angulated Channel JDEvolution Plus **		
EVA135C:	Tibase Interface Engaging H 10.0 for Angulated Channel JDEvolution Plus **		
EVA136C:	Tibase Interface Non Engaging H 10.0 for Angulated Channel JDEvolution Plus **		
EVCAITEC:	Interface CAD CAM Engaging for Conical Abutment JDEvolution Plus (H 7.5)	Torque recommended 15 Ncm	
EVA139C:	TiBase Engaging for Conical Abutment JDEvolution Plus (H 5.0)		
EVCAITC:	Interface CAD CAM Non Engaging for Conical Abutment JDEvolution Plus (H 7.3)		
EVCATIBC:	TiBase for Conical Abutment JDEvolution Plus (H 5.0)		
EVCAITAC:	Interface CAD CAM Non Engaging for Conical Abutment Angulated Channel (H 7.7)		
EVS110:	Prosthetic Screw Conical Abutment Direct	Torque recommended 15 Ncm	
EVS111:	Prosthetic Screw Conical Abutment Direct Long		
EVS112:	Prosthetic Screw for Angulated Channel Conical Abutment Direct ***		
EVANCN:	Implant Replica CAD CAM New JDEvolution Plus	Torque recommended 30 Ncm	
EVCAARCN:	Conical Abutment Replica CAD CAM New		
EVPMNC:	Pre-milled Abutment Blank JDEvolution Plus		

* for this screw use Angulated Screw Driver for JDTorque

** The EVS: prosthetic screw must be used for straight channels and the EVSA: prosthetic screw must be used for angled channels. Prosthetic screws are not included.


*** for screw use Angulated Screw Driver for Conical Abutment Direct JDTorque

Note: please contact info@jdentalcare.com for the digital library for Exocad, 3Shape and Dental Wings software


Conometric Abutment:

EVCP05C:	Conometric Abutment Straight H 0.5 JDEvolution Plus	Torque recommended 30 Ncm	
EVCP10C:	Conometric Abutment Straight H 1.0 JDEvolution Plus		
EVCP20C:	Conometric Abutment Straight H 2.0 JDEvolution Plus		
EVCP30C:	Conometric Abutment Straight H 3.0 JDEvolution Plus		
EVCP40C:	Conometric Abutment Straight H 4.0 JDEvolution Plus		
EV7CP20C:	7.5° Angulated Conometric Abutment H 2.0 JDEvolution Plus	Torque recommended 30 Ncm	
EVI5CP20C:	15° Angulated Conometric Abutment H 2.0 JDEvolution Plus		
EV22CP20C:	22.5° Angulated Conometric Abutment H 2.0 JDEvolution Plus		
EVA300C:	30° Angulated Conometric Abutment H 2.0 JDEvolution Plus	Torque recommended 30 Ncm	
EVA285C:	7.5° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA286C:	15° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA287C:	22.5° Angulated Conometric Abutment H 3.0 JDEvolution Plus		
EVA288C:	30° Angulated Conometric Abutment H 3.0 JDEvolution Plus	Torque recommended 30 Ncm	
EVCPS:	Screw for Conometric Abutment JDEvolution Plus		
EVA130C:	Cap With Pin for Intraoral Welding		
EVA131C:	Cap for Partial or Full-Arch Fixed Prosthesis		
EVA132C:	Anti-Rotation Cap for Single Element		
EVCPI C	Conometric Abutment Impression Coping PEEK		
EVCPAN	Conometric Abutment Replica		

Ball Abutment:

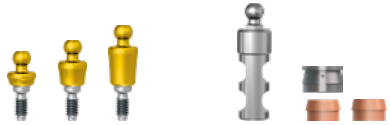EVBA15:	Ball Abutment H 1.5 JDEvolution Plus	Torque recommended 30 Ncm	
EVBA30:	Ball Abutment H 3.0 JDEvolution Plus		
EVBA50:	Ball Abutment H 5.0 JDEvolution Plus		
EVBA A	Ball Abutment Replica		
EVBAHC	Ball Abutment Elastic Retention Cap w/Container		
EVBA N	Ball Abutment Retention Cap		

Emi Abutment:

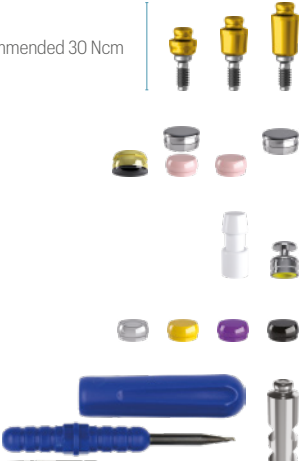EVEMI15:	Emi Abutment H 1.5 JDEvolution Plus	Torque recommended 30 Ncm	
EVEMI30:	Emi Abutment H 3.0 JDEvolution Plus		
EVEMI50:	Emi Abutment H 5.0 JDEvolution Plus		
EVEMISCK	Smart cap attachment Emi Abutment		
EVEMIH C	Emi Abutment Elastic Retention Cap w/ Container		
EVEMIN	Emi Abutment Elastic Retention Cap (1.2 kg)		
EVEMIH	Emi Abutment Container		
EVEMIIC	Emi Abutment Plastic Impression Coping		
EVEMIICS	Emi Abutment Steel Impression Coping		
EVEMINT	Emi Abutment White Retention Cap (1.8 kg)		
EVEMINY	Emi Abutment Yellow Retention Cap (0.6 kg)		
EVEMINP	Emi Abutment Purple Retention Cap (2.5 kg)		
EVEMINB	Emi Abutment Retention Cap for Laboratory		
EVEMIEIT	Extractor-insertion Tool For Caps		
EVEMIA	Emi Abutment Replica		

Implant replica:

EVAN: Implant Replica JDEvolution Plus



Bone Mill:

JDBMNC: Bone Mill Ø 5.0 and Guide JDEvolution Plus
JDBM6NC: Bone Mill Ø 6.0 and Guide JDEvolution Plus
JDBMGNN: Bone Mill Guide JDEvolution Plus



Short Drills:

JDDICS4: Short Drill Ø 4mm L 6 JDIcon Plus
JDDICS4D: Short Drill Ø 4mm L 6 Dense Bone JDIcon Plus
JDDICS5: Short Drill Ø 5mm L 6 JDIcon Plus



Implant and Prosthetic Drivers:

EVID: Implant Driver JDEvolution Plus
EVIDL: Implant Driver Long JDEvolution Plus
JDID127: Implant Depth Aligner JDEvolution Plus

EVSDPF15: Prosthetic Driver for JDTorque L 15 JDEvolution Plus
EVSDPF20: Prosthetic Driver for JDTorque L 20 JDEvolution Plus
EVSDPF25: Prosthetic Driver for JDTorque L 25 JDEvolution Plus
EVSDPF35: Prosthetic Driver for JDTorque L 35 JDEvolution Plus
EVSDPF60: Prosthetic Driver for JDTorque L 60 JDEvolution Plus


EVSDCAF: Conical Abutment Driver for JDTorque

EVSDP20: Prosthetic Driver L 20 JDEvolution Plus
EVSDP25: Prosthetic Driver L 25 JDEvolution Plus
EVSDP30: Prosthetic Driver L 30 JDEvolution Plus
JDPD105: Prosthetic Driver JDEvolution Plus for Surgical Driver

EVSDCA: Conical Abutment Driver

EVSDPF25A: Angulated Screw Driver for JDTorque L 25
EVSDPF30A: Angulated Screw Driver for JDTorque L 30
EVSDPF25AV: Angulated Screw Driver for Conical Abutment Direct JDTorque L 25
EVSDPF30AV: Angulated Screw Driver for Conical Abutment Direct JDTorque L 30
EVSDP32AV: Screw Driver for Conical Abutment Direct EVS112 (Straight Channel)

EVSUD: Surgical Driver JDEvolution
JDG: JDGuide
EVGMC: Implant Mounter for JDEvolution
EVGMLC: Implant Mounter Long JDEvolution
JDID102: ZAPA Tool
JDSRP: Dental Implant Depth Probe

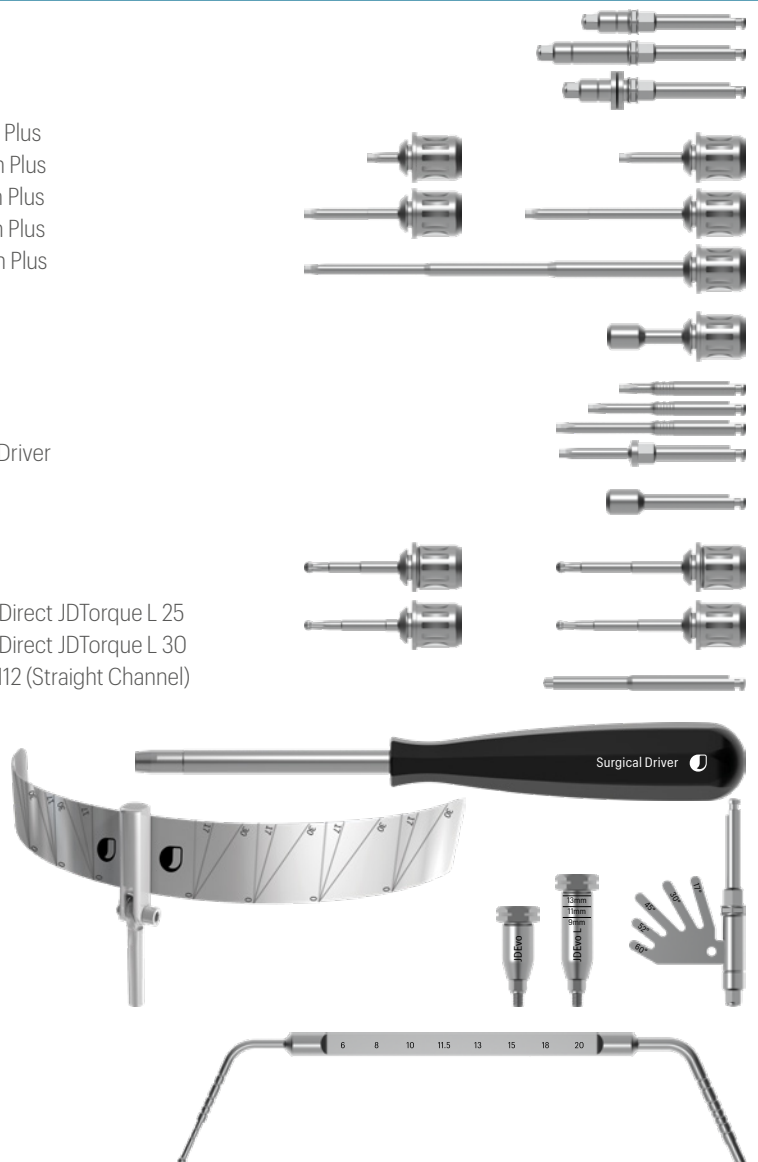

Direction indicators:

JDDI: Direction Indicator
JDDIS: Direction Indicator Short
JDDI17: Direction Indicator 17°
JDDI30: Direction Indicator 30°
JDDI45: Direction Indicator 45°



DRILLING PROTOCOLS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants

	HEALED BONE		POST EXTRACTIVE BONE		 Short Implant L 6mm JDEvolution Plus insertion in maxilla	
IMPLANT DIAMETER	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE		
site preparation in maxilla						
Ø 3,7	2,0 2,4 2,8 up to the 1 st laser mark L6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L6mm	2,0 2,4 2,8 at the entrance	2,0 2,4 2,8 at the entrance	 Short Implant L 6mm JDEvolution Plus insertion in mandible	
Ø 4,0	2,0 2,4 2,8 up to the 1 st laser mark L6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L6mm	2,0 2,4 2,8 at the entrance	2,0 2,4 2,8 at the entrance		
Ø 4,3	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 up to the 1 st laser mark L6-8mm 3,6 at the entrance	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 at the entrance		
Ø 5,0	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm		
Ø 6,0	2,0 2,4 2,8 3,2 3,6	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0	2,0 2,4 2,8 3,2 3,6 4,0		
site preparation in mandible					 Short Implant L 6mm JDEvolution Plus insertion in maxilla	
Ø 3,7	2,0 2,4 2,8 3,2 up to the 2 nd laser mark L 8mm 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 at the entrance	2,0 2,4 2,8 3,2 at the entrance		
Ø 4,0	2,0 2,4 2,8 3,2 up to the 2 nd laser mark L 8mm 3,6 up to the 1 st laser mark L 6mm 4,0 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6mm 3,6 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 up to the 1 st laser mark L 6mm 3,6 up to the 1 st laser mark L 6mm		
Ø 4,3	2,0 2,4 2,8 3,2 3,6 up to the 2 nd laser mark L 8mm 4,0 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 up to the 2 nd laser mark L8mm 4,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 at the entrance	2,0 2,4 2,8 3,2 3,6 at the entrance		
Ø 5,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 up to the 2 nd laser mark L 8mm 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance	2,0 2,4 2,8 3,2 3,6 4,0 at the entrance		
Ø 6,0	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 1 st laser mark L 6mm 5,4 up to the 1 st laser mark L 6mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4 4,8 up to the 2 nd laser mark L 8mm 5,4 up to the 2 nd laser mark L 8mm	2,0 2,4 2,8 3,2 3,6 4,0 4,4	2,0 2,4 2,8 3,2 3,6 4,0 4,4		


Note: All measurements in mm

DRILLING PROTOCOLS - Ø 4,0 T IMPLANTS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants.

IMPLANT DIAMETER	HEALED BONE		POST EXTRACTIVE BONE		
	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE	
site preparation in maxilla					Short Implant L 6mm JDEvolution Plus insertion in maxilla
Ø 4,0	2.0 2.4 2.8	2.0 2.4 2.8 3.2 up to the 2 nd laser mark L 8mm	2.0 2.4 2.8 up to the 2 nd laser mark L 8mm	2.0 2.4 2.8	short Ø 4,0 L 6 - Start the osteotomy with standard twist drill Ø 2.0mm and Ø 2.4mm - Complete with the Ø 4mm L 6mm drill code JDDICS4
site preparation in mandible					Short Implant L 6mm JDEvolution Plus insertion in mandible
Ø 4,0	2.0 2.4 2.8 3.2 up to the 2 nd laser mark L 8mm 3.6 up to the 2 nd laser mark L 8mm	2.0 2.4 2.8 3.2 3.6 4.0 up to the 1 st laser mark L 6mm	2.0 2.4 2.8 3.2 up to the 2 nd laser mark L 8mm 3.6 up to the 2 nd laser mark L 8mm	2.0 2.4 2.8 3.2 3.6 up to the 2 nd laser mark L 8mm	short Ø 4,0 L 6 - Start the osteotomy with standard twist drill Ø 2.0mm, Ø 2.4mm and Ø 2.8mm - Complete with the Ø 4mm L 6mm drill code JDDICS4D

Note: all measurements in mm

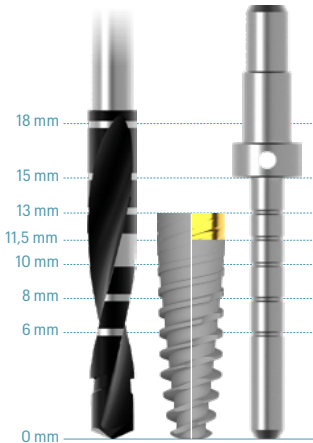



JDENTALCARE

www.jdentalcare.com

JDentalCare, the JDentalCare logo and other trademarks mentioned in this document are, excluding specifically mentioned exceptions, trademarks property of JDentalCare.

The images of the products represented in this brochure are not necessarily to scale.

Note: Not all products manufactured by JDentalCare and represented in this brochure have been regulatory cleared/released for sale in all countries and/or have been licensed in accordance with FDA Rules and Regulations.

For more information about current product assortment and availability, please contact your local JDentalCare distributor. For prescription only.

Caution: Federal (US) law restricts this device to sale by or on the other for a licensed dentist.

See instructions for use for more information such as indications, warnings, precaution and contraindications.

EVPCAT
REV.71 of 2026-09-11